


FITTING INSTRUCTIONS: Thread top strap through buckles. Repeat for bottom strap and buckle.




1. Position bottom strap around the back of the head, just below the ears.



2. Position top strap around the back of the head just above the ears.

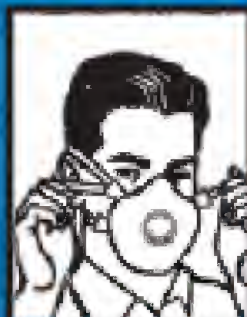

3. Adjust the tension by pulling on end tab of each strap.

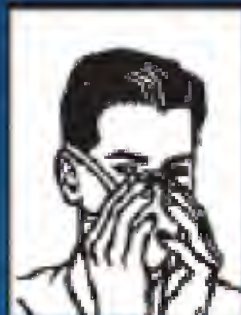

4. Using both hands, mold the metal nosepiece to the shape of your nose.



5. To release strap tension without removing the respirator, push out on back of buckle.



6. Check seal of respirator on the face, prior to use; for proper fit. A) Cup both hands over the front of the respirator being careful not to disturb position. B) Inhale sharply. A negative pressure should be felt inside respirator. If leakage is suspected, adjust position of respirator and/or strap and repeat procedure.

NOTE: If a proper fit cannot be achieved, do not enter contaminated area. See your supervisor.